


## Charging Light

- RED when charging
- GREEN when charged



## Unit Functions

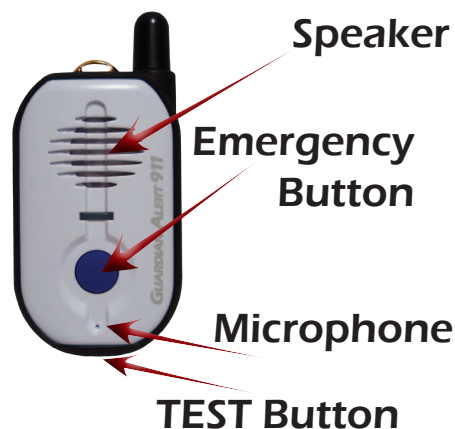

## Charging the Pendant



- Fully charge pendant by putting it in the charging cradle prior to use. It may take 3 to 4 hours to ensure a full charge.
- The Charging Light on the pendant will be RED when charging and turn GREEN when fully charged. We recommend that the pendant be charged monthly.

## Checking your System

Push the **WHITE** test button on the bottom of the pendant for **2 seconds** to **test the battery**. Then, press it **again** for **4 seconds** to **test the system** and listen for voice prompt from the pendant:



- "All Systems are OK" is heard if the system is correctly installed.
- If any other voice prompt is heard, the unit will describe the problem:
  - a. "Cannot detect telephone line" - Plug telephone cord into a working phone jack
  - b. "Cannot detect dial tone" - An extension phone is off the hook.
  - c. "Recharge battery soon" - Place pendant in charging base for 3-4 hours.
  - d. "Base and pendant out of range" - There is no power to the base, or the base and pendant are not communicating.

**We recommend that you perform this System Test weekly.**